

Qualitix® sterilization pouches

consumables

Fitting any small to mid-size laboratory, medical, dental or veterinary device, the Qualitix® single use pouches provide easy and effective wrapping of items requiring sterilization. Made with robust medical grade paper, the pouches tightly close with self-sealing strip to maintain sterility. Content remains always visible through PET / CPP transparent laminated film.



Ask for free samples


The Qualitix® sterilization pouches

Product advantages

- Robust durable material
- Quick content identification
- Large self-sealing strip, easy to close
- Chevron top seal for easy opening
- Suitable for steam and EO gas sterilization
- Two sterilization indicators
- Five different sizes
- Conformity to ISO 17665-06, ISO 11135, CE




Self-sealing - quick opening ①
Simply peel off protection strip, fold and press with fingers. Pouch remains tightly sealed to maintain product sterility.

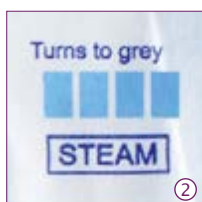

Sterilization indicators ②
Two chemical indicators included to confirm proper sterilization process. Blue indicator turns dark grey after autoclaving, red indicator turns to yellow after ethylene oxide (EO) gas sterilization.



Integrity control ③
Blue stained laminated film enables easy control of film/pouch integrity. Do not use a damaged pouch.

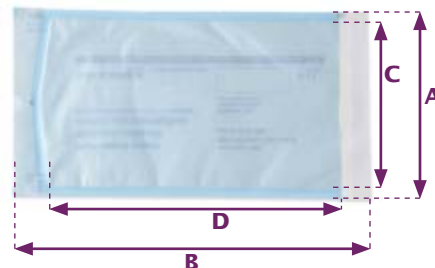

Ordering information

Description	Dimensions		Material	Packaging	Cat. No.
Pouch, self-sealing	83 x 160 mm	3.3" x 6.3"	Med. grade paper/PET/CPP	200 / pk	3210.0816
Pouch, self-sealing	90 x 260 mm	3.5" x 10.2"	Med. grade paper/PET/CPP	200 / pk	3210.0926
Pouch, self-sealing	135 x 260 mm	5.3" x 10.2"	Med. grade paper/PET/CPP	200 / pk	3210.1326
Pouch, self-sealing	190 x 360 mm	7.5" x 14.2"	Med. grade paper/PET/CPP	200 / pk	3210.1936
Pouch, self-sealing	300 x 395 mm	11.8" x 15.6"	Med. grade paper/PET/CPP	200 / pk	3210.3040

Dimensions

A (mm)	B (mm)	C (mm)	D (mm)
83	160	65	110
90	260	70	205
135	260	115	205
190	360	170	300
300	395	270	335